AbViser® AutoValve™

Intra-Abdominal Pressure Monitoring Device

Manufactured by:

Wolfe Tory Medical, Inc.

79 West 4500 South, Suite 16
Salt Lake City, Utah 84107
Tel: 801-281-3000 Fax: 801-281-0708
Toll Free: 1-888-380-9808
www.wolfetory.com

EU Authorized Representative:

Cavendish Scott Ltd.
P.O. Box 107
SG5 1 FW, England

Elegant solutions for complex issues.

AbViser is a registered trademark and AutoValve is a trademark of Wolfe Tory Medical, Inc.
TruWave is a registered trademark of Edwards Lifesciences Corporation.
Transpac is a registered trademark of Hospira, Inc.

Made in the U.S.A. US and Foreign Patents Pending

P/N 30-0212, REV. 06/06

AbViser® AutoValve™

Intra-Abdominal Pressure Monitoring Device

Instructions For Use

REF

Catalog Numbers:

ABV300, ABV301

ABV320, ABV321

ABV330, ABV331

WOLFE TORY
MEDICAL, INC.

WOLFE TORY
MEDICAL, INC.

AbViser® AutoValve™

Intra-Abdominal Pressure Monitoring Device

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Bruksanvisning

REF ABV300, ABV301, ABV320, ABV321, ABV330, ABV331

AbViser® AutoValve™ Intra-Abdominal Pressure Monitoring Device

Indications For Use:

The AbViser AutoValve device is indicated for:

- Infusing fluid into a urinary bladder catheter.
- Measuring the hydrostatic pressure within the urinary bladder.

IMPORTANT:

- Models ABV300 and ABV301 of the AbViser AutoValve Intra-Abdominal Pressure Monitoring Device do **NOT** include a pressure transducer. They will adapt to the pressure transducer and patient monitor used in your ICU. Models ABV320, ABV321, ABV330 and ABV331 **DO** include a pressure transducer.
- The AbViser Intra-Abdominal Pressure Monitoring Device may be used on a patient for as long as the urinary catheter and urinary drain bag are in place. The AbViser must be replaced at the time the urinary catheter and/or urinary drain tubing is replaced.

CAUTION:

Luer lock connections can come loose during shipping and handling. Be sure to check that all luer lock connections are tight prior to use.

Single Patient Use Only


Rx Only: Federal (U.S.A) law restricts this device to sale by or on the order of a physician.


Latex Free


Sterile unless package is damaged or open


AbViser AutoValve Sterile Device Package Contents:

- AbViser AutoValve
- 20 mL 3-ring syringe with protective silicone sleeve
- Double check valve
- Tubing and spike for sterile saline bag connection
- Infusion and pressure transducer tubing
- Zeroing stopcock
- Edwards Lifesciences TruWave® Transducer (Models ABV320 & ABV321)
- Hospira Transpac® IV Transducer (Models ABV330 & ABV331)
- Caps for dead-ending transducer (Models ABV300 & ABV301)
- Tape - to provide extra security between urinary catheter and AbViser valve
- S-hook
- Sterile drape
- Instructions for use

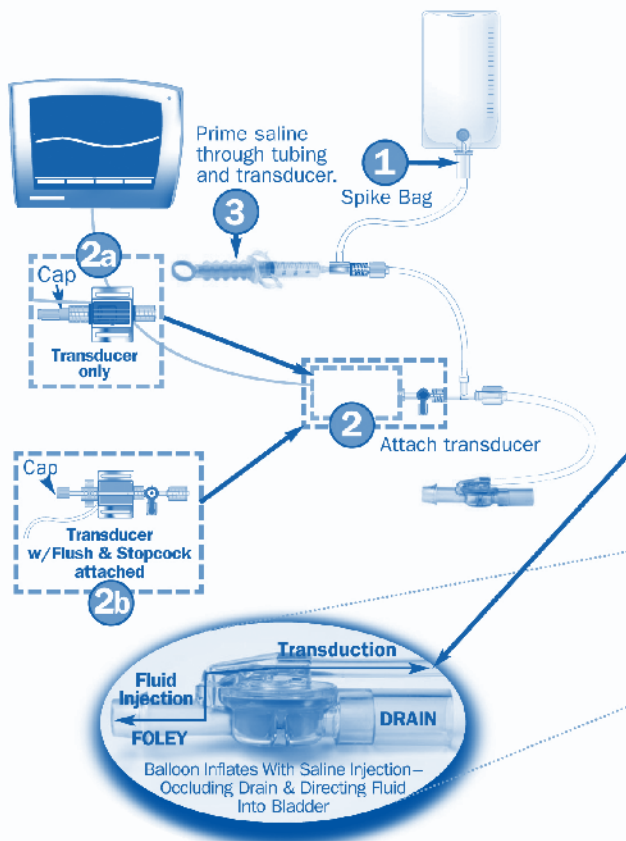AbViser[®] AutoValve[™] Intra-Abdominal Pressure Monitoring Device

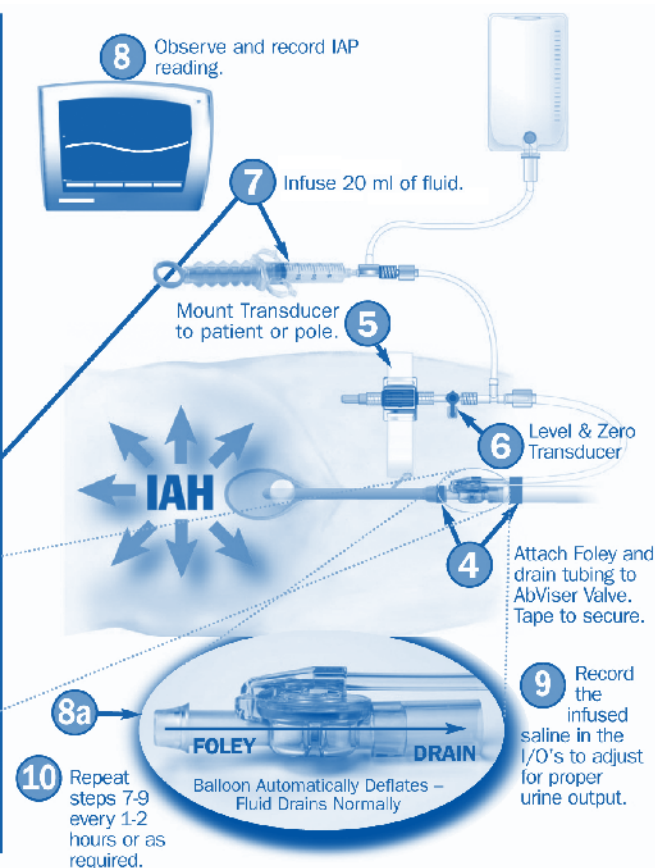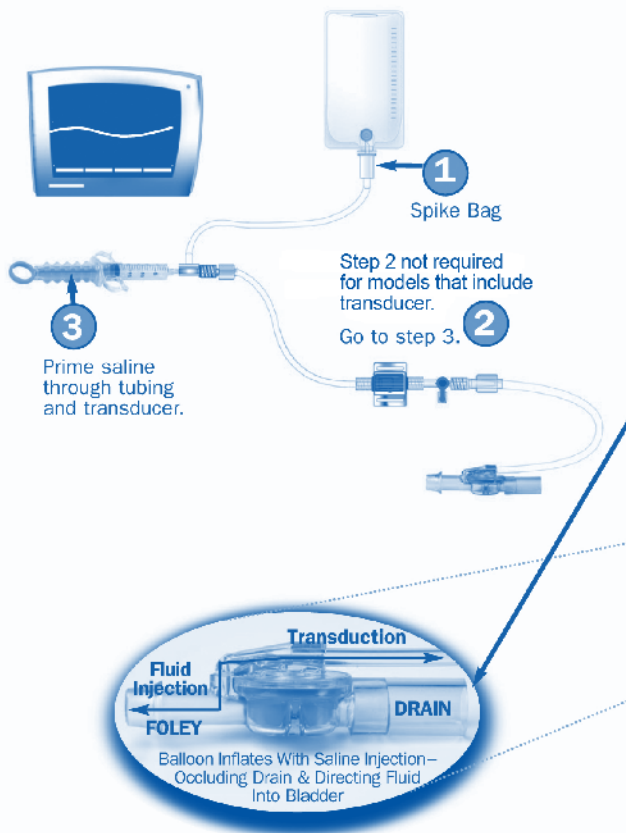AbViser® AutoValve™ Intra-Abdominal Pressure Monitoring Device

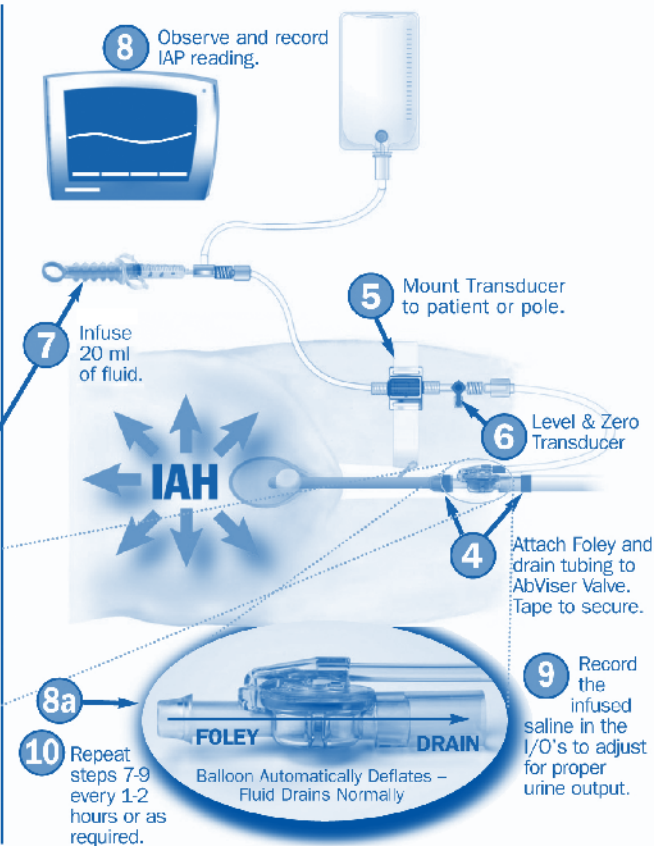AbViser® AutoValve™ Intra-Abdominal Pressure Monitoring Device

Figure 1 - ABV300, ABV301

Figure 2 - ABV320, ABV321, ABV330, ABV331

1 Assembling the AbViser

- 1.1 Hang a bag of sterile saline on an IV pole. A PRESSURE CUFF IS NOT REQUIRED.
- 1.2 Prep the urinary catheter/drain tubing connection with antiseptic solution (you will be disconnecting them and inserting the AbViser valve between them).
- 1.3 Open the AbViser tray.
- 1.4 Grasp the saline bag tubing spike. Remove protective cap and insert the spike into the saline bag. For AbViser AutoValve models ABV300 and ABV301 you will need to open a transducer as described below, for models ABV320, ABV321, ABV330 and ABV331 the transducer is included and you can proceed directly to Step 3.

2 Attaching the Transducer (Models ABV300 and ABV301):

- 2.1 Open a pressure transducer kit (not included) in sterile fashion and drop contents onto the open AbViser tray.
NOTE: The AbViser adapts to the pressure transducer used in your ICU. Many pressure transducer variations exist, all which are compatible with the AbViser. Although a flush device is not required, the AbViser can be used with a flush device. OPTION 2 below is available because many flush devices are permanently attached or difficult to disconnect from the transducer.
- 2.2 OPTION 1: Using sterile techniques, remove all attachments from the transducer (Flush Device, Drip Chamber, Tubing, Stopcocks, Etc) until the transducer has only luer connections. OPTION 2: Remove attachments per OPTION 1 except leave flush device and/or stopcock attached. Be sure the flush device is attached distally and capped at the end. (See Figure 1)
- 2.3 Connect the pressure transducer to the stopcock on the AbViser tubing assembly. Cap the opposite end of the transducer (Figure 1) or flush device.

** IMPORTANT NOTE **

The AbViser DOES NOT require a pressure cuff around the saline bag.

4 Priming the AbViser

- 3.1 Aspirate and compress the syringe to flush all the air from the infusion tubing until saline runs out of the AbViser AutoValve and all the bubbles are out of the tubing.
- 3.2 For models ABV300 and ABV301 only, remove the end cap from the transducer, open the stopcock and compress the syringe until saline runs out the end of the transducer, then recap the transducer. Other models do not require this step.

3 Connecting the AbViser to the Patient

- 4.1 Place the sterile drape under the patient's urinary catheter/urinary drain bag connection.

- 4.2 Prep the connection of the urinary catheter and urinary drain bag with antiseptic solution (if not already done in Step 1).
- 4.3 Detach the Foley catheter from the drain connection.
- 4.4 Pick up the AbViser valve. Tear perforation of protective bag, exposing barbed connector. Push the urinary catheter over the barb. Slide bag along tubing exposing the urinary drain connector. Connect drain tubing to AbViser valve. Tear plastic bag at second perforation and discard.
- 4.5 Ensure catheter connection is dry, and then wrap a strip of tape (included in kit) around the end of the catheter to keep it snugly attached to the AbViser valve. (Figures 1 & 2) Wrap a strip of tape (included in the kit) around the AbViser valve infusion tubing and the urinary drain tubing – this will prevent kinks from occurring in the infusion tubing. **CAUTION: Once the AbViser valve is connected to the urinary drain tubing, it cannot be disconnected without damaging the valve.**
- 4.6 Hang the syringe at a convenient spot using the S-hook.

5 Attach the pressure transducer to the patient or the pole.

- 5.1. Connect the pressure transducer cable to the pressure transducer.

6 Zero the transducer to atmosphere at the level of the mid-axillary line. After zeroing the transducer, ensure the stopcock is open to the transducer.

7 Measuring Bladder Pressure with the AbViser

- 7.1 Patient must be supine (flat on back).
- 7.2 Retract the plunger until 20 mL of fluid is within the syringe. Compress the syringe plunger over several seconds, infusing the fluid into the bladder.

8 Allow the system to equilibrate and then note the pressure reading on the monitor at end expiration. This IAP reading will last approximately 1 minute, at which point the valve will automatically open (drain) and the pressure will rapidly decrease to zero. (Note: If the pressure does not completely decrease to zero, there may be fluid remaining in the drainage tubing. Drain the fluid into the drainage bag and the reading should go to zero, indicating the AbViser AutoValve has opened into the drain position.) The valve will remain open (in drain position) until the next volume of fluid is infused.

9 Record the volume of fluid infused in the fluid I&O's sheet.

10 Repeat step 7-9 every 1-2 hours or as directed by your hospital protocol.

AbViser® AutoValve™ Intra-Abdominal Pressure Monitoring Device

Wolfe Tory Medical, Inc.
79 West 4500 South, Suite 16
Salt Lake City, Utah 84107
Tel: 801-281-3000 Fax: 801-281-0708
Toll Free: 1-888-380-9808 www.wolfetory.com

Elegant solutions for complex issues.

Made in the U.S.A.
US and Foreign Patents Pending

AbViser is a registered trademark and AutoValve is a trademark of Wolfe Tory Medical, Inc.

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Tillverkas i U.S.A.
AbViser är ett registrerat varumärke och AutoValve är ett varumärke som tillhör Wolfe Tory Medical, Inc.